


NIPPY S+

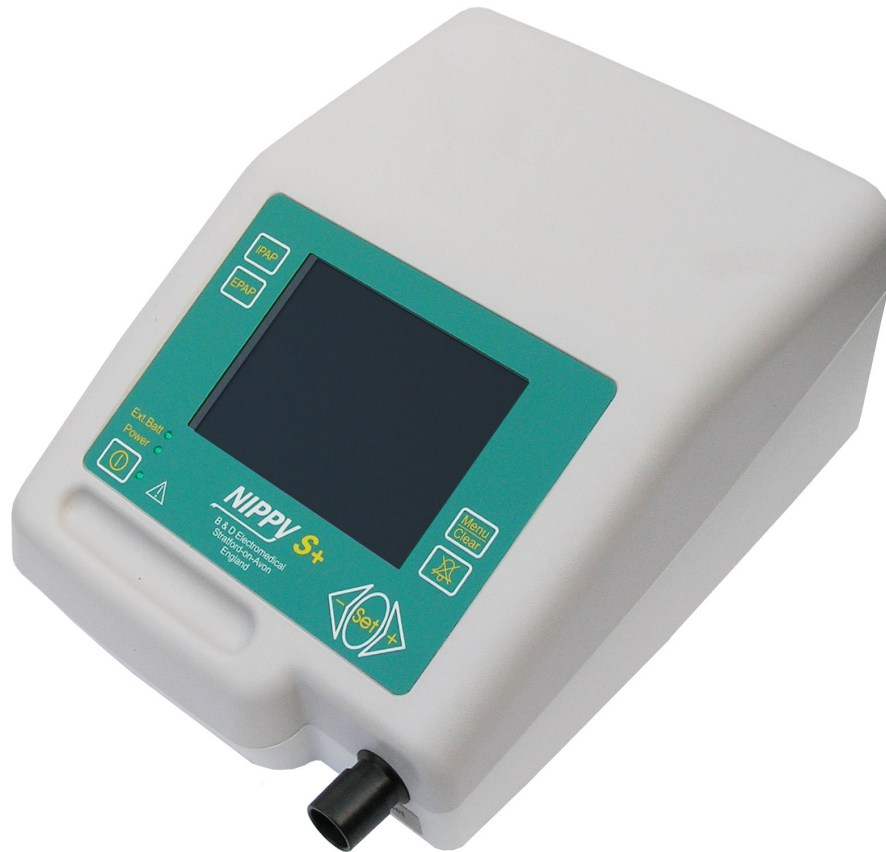

Bi-Level positive pressure ventilator

for use with Nasal Mask or Nasal Pillows

Optional INTERNAL BATTERY

Run times up to 12 hours

The **Nippy S+** is a pressure controlled positive pressure ventilator employing state of the art electronics, to give consistent and accurate support for patients requiring non invasive ventilation via a nasal mask or nasal pillows.

The large colour LCD screen gives unrivalled clarity when setting up or reading stored data.

Intuitive software allows for easy setting up, with on-screen help available for each function.

Pressure Support Ventilation

Settings and airway measurements are displayed on-screen with more advanced features accessed via an easy to follow menu.

Compliance data may be viewed on-screen.

Full event log may be downloaded to a PC.

The **Nippy S+** is equipped with a universal mains power inlet, a dc power inlet and an

Tel: 01789 293460 Sales.
Fax: 01789 262470


01789 721577 Technical Support
www.nippyventilator.com

internal battery. Weighing only 3.6kg,
world-wide travel has never been easier.

Tel: 01789 293460 Sales.
Fax: 01789 262470

01789 721577 Technical Support
www.nippyventilator.com

Nippy S+

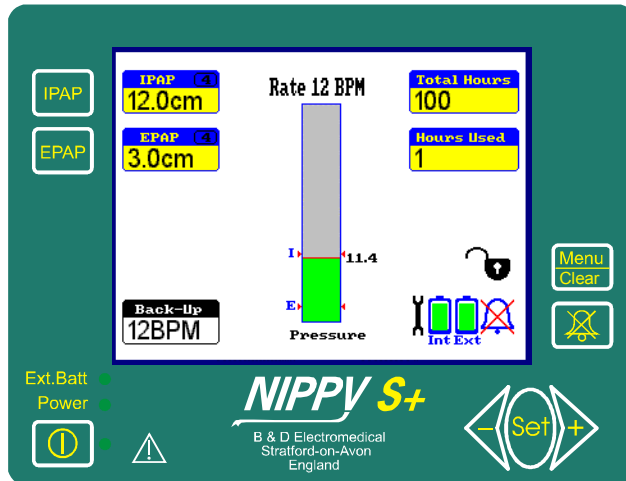

The Nippy S+ is intended for use by patients who require respiratory support via nasal mask or nasal pillows.

For the treatment of:-
Chronic Obstructive Pulmonary Disease (COPD)

Decreased drive to breathe (Central Hypoventilation)

The Nippy S+ is equally suited to hospital or domicilliary use.

Technical Specifications

Supply Voltage	100 - 240 Vac
Dimensions cm (LxWxH)	30 x 22 x 13
Weight	3.6 kg
With battery	4.5kg
Pressure Range IPAP	3 – 38 cm H ₂ O
Pressure Range EPAP	3 – 20 cm H ₂ O
Fixed Back-up Rate	12BPM
Ti	1.5 sec
Optional Internal battery	24Vdc 6.8 Ahr
Running time	4 -12 hours
Optional External battery	24Vdc 6.8 Ahr
Running time	4 -12 hours
Alarms	Mask off, Power Failure

Features

- ◆ Pressure Support – Spontaneous breathing augmented by the ventilator
- ◆ User friendly intuitive software
- ◆ Large colour LCD display, clearly shows all settings.
- ◆ Lockable settings.
- ◆ Alarms with easily understood on-screen messages.
- ◆ On-screen compliance data.
- ◆ Full event log showing all adjustments and alarm events may be downloaded to a PC.
- ◆ Universal mains input and internal battery operates anywhere in the world.
- ◆ Adjustable fast response flow triggers.
- ◆ Can be connected to Nurse call systems using the optional nurse call adapter.

Tel: 01789 293460 Sales.
Fax: 01789 262470

01789 721577 Technical Support
www.nippyventilator.com